

Hydrofluoric Acid Content of Lithium-ion Battery Electrolyte

Potentiometric Titration in the Mettler Toledo Titration Excellence T5 To Determine HF Content

Lithium-ion batteries (LIB) are commonly used for portable electronic devices and electric vehicles. They are growing in popularity because of their high energy density. The lithium-ion battery stores energy that is released by an electrochemical reaction between the anode and the cathode material. Both the anode and the cathode contain positively charged lithium ions. During discharge, the oxidation reaction at the anode produces electrons and lithium cations. The electrons flow through an external wire to the cathode. To complete the electric circuit, lithium cations, flow through the electrolyte to the cathode where they recombine with the electron in a reduction reaction. Although the lithium ions formally do not partake in the electrochemical reaction, it is important that they are efficiently carried to the respective electrode by the electrolyte.

Keywords or phrases: Lithium-ion, electrolyte, hydrofluoric acid, titration, battery, lithium, Mettler Toledo

Introduction

The electrolyte is usually composed of the electrolyte salt and a number of additives (e.g. organic carbonates). Lithium hexafluorophosphate LiPF_6 is the most frequently used electrolyte salt in LIB electrolytes. It dissolves well in organic carbonates, is electrochemically stable and helps to conduct lithium cations. However, the use of LiPF_6 also brings up a problem. It is not stable against hydrolysis. Traces of water in the electrolyte lead to the formation of hydrofluoric acid HF, which is considered as battery poison destroying the electrodes and lowering the capacity. Therefore, the amount of HF in the

Application Showcase

lithium-ion battery electrolyte should be as low as possible. The HF content serves as a quality indicator and is listed in the specification of the electrolyte.

This application note describes a method to determine the hydrofluoric acid content in a lithium-ion battery electrolyte by an aqueous acid-base titration with potentiometric indication.


Sample Preparation and Procedures

Add approximately 40 mL of deionized water into a 100 mL titration beaker and place the beaker on the balance.

Tare the balance and add ca. 3 g of the electrolyte using a syringe.

Then place the beaker on the InMotion autosampler rack and start method as shown in Table 2.

Chemistry

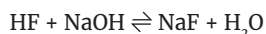

Chemicals

Titrant: Aq. sodium hydroxide, NaOH, $c(\text{NaOH}) = 0.01$ mol/L.

Sample: Lithium hexafluorophosphate solution in ethylene carbonate and ethyl methyl carbonate, 1.0 M, LiPF_6 in EC/EMC=50/50 (v/v), battery grade

Analyte: Hydrofluoric acid, HF, $M = 20.01$ g/mol, $z = 1$.
Auxiliary Reagents: deionized water.

Standard: Potassium hydrogen phthalate, $\text{C}_8\text{H}_5\text{KO}_4$,
 $M = 204.22$ g/mol, $z = 1$

Instruments and Accessories

- › Titration Excellence T5
- › InMotion flex 100 mL
- › Burette DV1010 10 mL
- › Compact Stirrer
- › DGi111-SC or DGi115-SC electrode

› Excellence balance XPR603S Delta Range

Results

	Sample Size (g)	R1: Content (ppm)
1	2.8542	75.02
2	3.1214	76.78
3	2.897	78.62
4	2.7795	75.61
5	2.8144	76.52
6	3.2041	76.97
Mean		76.59
s		1.24
srel		1.62%

The results show a good repeatability with a relative standard deviation of 1.7%. The results are listed in chronological order. The inexistence of a trend towards increasing values indicates that the hydrolysis during the sample handling is negligible.

Measured Values

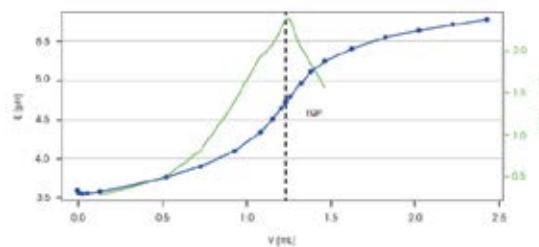


Figure 1: Titration curve (blue) and the first derivative of the titration curve (green) of sample 6/6.

Application Showcase

Method

Sample		Evaluation and Recognition	
Sample Type	Sample	Procedure	Standard
Number of IDs	1	Threshold	1.5 pH/mL
ID 1	Electrolyte	Tendency	Positive
Entry type	Weight	Ranges	0
Lower Limit	0	Add. EQP criteria	No
Upper Limit	5 g	Termination	
Density	1 g/mL	At Vmax	5 mL
Correction Factor	1	At potential	No
Temperature	25.0 °C	At slope	No
Titration reader	None	After number of recognised EQPs	Yes
Number of sample factors	0	Number of EQPs	1
Titration stand		Combined Termination Criteria	No
Type	InMotion T/Tower A	Condition	No
Titration Stand	InMotion T/1A	Calculation R1	
Head Position	Sample	Result	Content
Stir		Result Unit	ppm
Speed	30%	Formula	$R1 = Q \cdot C / m$
Duration	5 s	Constant C	$C = M \cdot 1000 / z$
Titration (EQP)1		M	M[Hydrofluoric acid]
Titration		Decimal Places	3
Titrant	NaOH	Result Limits	No
Concentration	0.01 mol/L	Extra Statistical Functions	No
Sensor		Send to Buffer	No
Type	pH	Write to Smart Tag	None
Sensor	DGi111-SC	Condition	No
Unit	mV	Rinse	
Temperature Acquisition		Auxiliary Reagent	Water
Temperature Measurement	No	Rinse Cycles	1
Stir		Vol. per cycle	10 mL
Speed	30%	Position	Current Position
Predispense		Drain	No
Mode	None	Condition	No
Wait time	0 s	End of Sample	
Control		Park	
Control	Cautious	Titration Stand	InMotion T/1A
Mode	Acid/Base	Position	Condition Beaker
Titrant Addition	Dynamic	Condition	No
dE (set value)	10 mV	Record	
dV (min)	0.005 mL	Report Template	Titration Report
dV (max)	0.2 mL		
Mode	Equilibrium controlled		
dE (set value)	0.5 mV		
dt	2 s		
t (min)	3 s		
t (max)	30 s		

Table 2: Experimental Method

Application Showcase

t (s)	V (mL)	E (pH)	dE/dV (mV/mL)
0	0	3.586	n/a
3	0.005	3.553	n/a
6	0.01	3.55	n/a
9	0.0225	3.549	n/a
12	0.0535	3.554	n/a
17	0.131	3.58	0.3
23	0.325	3.66	0.38
30	0.525	3.766	0.51
37	0.725	390.10%	0.82
47	0.925	4.097	1.38
60	1.083	4.343	1.93
72	1.153	4.505	2.07
86	1.204	4.646	2.25
EQP	1.23252	4.726	2.37
100	1.2555	4.7	2.37
118	1.318	4.964	2.05
138	1.3805	5.113	1.86
159	1.4645	5.25	1.86
181	1.623	540.80%	n/a
206	1.823	555.50%	n/a
222	2.023	5.643	n/a
238	2.223	5.723	n/a
251	2.423	5.783	n/a

Conclusion

LiPF_6 hydrolyses slowly in aqueous solutions, while it is hydrolysing faster in organic solutions with traces of water¹. Thus, the titration is conveniently performed in water using aqueous NaOH as titrant.

To prevent the further formation hydrofluoric acid by hydrolysis, the sample should be titrated readily.

This method was only tested on the presented sample. The compositions of Lithium-ion battery electrolytes vary, and some electrolytes may contain additives that hydrolyse to form of hydrofluoric acid very quickly. For such samples, it may be necessary to strictly work in non-aqueous conditions.

The titre of aq. 0.01 M sodium hydroxide solution is determined using potassium hydrogen phthalate (KHP) as the primary standard. As the optimum amount of standard per measurement is very low (5 mg), a

5 mg/mL aq. KHP standard solution is prepared.

After every sample, the electrode, stirrer and tubing are rinsed with deionized water.

The titrations are performed in plastic beakers as HF is able to etch glass.

Waste Disposal and Safety Measures

Although the hydrofluoric acid content of the electrolyte is quite low, every contact with skin or eyes should be prevented. Use gloves, safety goggles and a lab coat. Work in a fume hood. After the titration, all solutions are neutralized and disposed as aqueous waste.

Reference

1. Stich M., Göttlinger M., Kurniawan M., Schmidt U., Bund A., J. Phys. Chem., 2018, 122, 16, 8836–8842.

Titration Excellence T5

